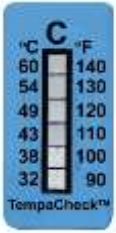
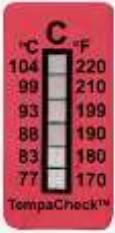


## Instructions for Use: TempaChek™ Irreversible Thermometers

<b>Brand Name of Product</b>	TempaChek™-90 & TempaChek™-170 Irreversible Thermometers
<b>Generic Name of Product</b>	Irreversible thermometers
<b>Product Code Number(s)</b>	Temp-90, Temp-170
<b>Purpose of Product</b>	To monitor key temperatures during the wash cycle of an automated instrument washer. <ul style="list-style-type: none"> <li>• Temp-90 – monitoring of the cold water (pre-wash) rinse.</li> <li>• Temp-170 – monitoring of the thermal disinfection rinse.</li> </ul>
<b>Range of Applications for Product</b>	Automated instrument washers
<b>Key Specifications of Product</b>	<ul style="list-style-type: none"> <li>• Irreversible thermometer – records the highest temperature achieved during exposure.</li> <li>• No residue adhesive allows placement on a clean, dry, smooth flat surface.</li> <li>• Easy to interpret color change from silver to black.</li> <li>• Graduated temperature scale to record range of temperatures.</li> <li>• Single-use.</li> <li>• Precision to +/-1°C.</li> </ul>

<b>Shipping &amp; Storage</b>	
<b>Shipping Conditions &amp; Requirements</b>	Temp-90 – do not allow ship temperature to exceed 85°F. Temp-170 – do not allow ship temperature to exceed 165°F.
<b>Storage Conditions</b>	Keep in cool, dry place out of direct sunlight. Temp-90 – do not allow temperature to exceed 85°F. Temp-170 – do not allow temperature to exceed 165°F.
<b>Packaging Contents</b>	10 per pack
<b>Shelf Life</b>	Expires two years from the date of manufacturer.

<b>Instructions for Using Product</b>	
<b>Description of Use(s)</b>	For monitoring key temperatures during the wash cycle of an automated instrument washer.
<b>Preparation for Use</b>	Be sure surface for adherence is dry, clean and smooth.
<b>Diagrams (drawings, pictures)</b>	<div style="display: flex; justify-content: space-around; align-items: center;"> <div style="text-align: center;">  <p>TempaChek™-90</p> </div> <div style="text-align: center;">  <p>TempaChek™-170</p> </div> </div>
<b>Steps for Use of Product</b>	<ul style="list-style-type: none"> <li>• TempaChek™-90             <ol style="list-style-type: none"> <li>1. Use one (1) TempaChek™-90 per washer.</li> <li>2. Peel thermometer from release paper.</li> <li>3. Apply to any clean, dry flat surface ensuring that the indicator has adhered to the surface.</li> <li>4. TempaChek™-90 should be read during the cold water rinse (pre-wash) stage.</li> <li>5. Record results on log sheet.</li> <li>6. Temperature should not exceed 100°F. Immediately report any result that exceeds this temperature.</li> </ol> </li> <li>• TempaChek™-170             <ol style="list-style-type: none"> <li>1. Use one (1) TempaChek™-170 on each level of the instrument rack.</li> <li>2. Peel thermometer from release paper.</li> <li>3. Apply to any clean, dry flat surface ensuring that the indicator has adhered to the surface.</li> <li>4. TempaChek™-170 should be read after the thermal disinfection stage and before the drying stage.</li> <li>5. Record results on log sheet.</li> </ol> </li> </ul>

	6. Report any deviation from targeted temperature.
<b>Interpretation of Results</b>	
<b>Contraindications of Test Results</b>	
<b>Documentation</b>	
<b>Special Warnings and Cautions</b>	<ol style="list-style-type: none"> <li>1. Peel off the thermometer(s) at the end of the wash cycle.</li> <li>2. Surfaces may be hot, so take caution and wear appropriate safety gear.</li> <li>3. Allowing the thermometer to stay on the surface for an extended period and for repeated cycles may make it very difficult to remove and adhesive residue may remain.</li> <li>4. If adhesive residue remains on the surface, may be removed with acetone or isopropyl alcohol solvent.</li> <li>5. Apply TempaChek™-90 and TempaChek™-170 labels only to a flat surface. Do not apply product to a round surface (the wash rack) as it may increase the risk of peeling and falling off during the wash cycle.</li> </ol>
<b>Disposal</b>	May be disposed of in regular trash.

<b>Reprocessing Instructions</b>	
<b>Point of Use</b>	
<b>Preparation for Decontamination</b>	
<b>Disassembly Instructions</b>	
<b>Cleaning – Manual</b>	
<b>Cleaning – Automated</b>	
<b>Disinfection</b>	
<b>Drying</b>	
<b>Maintenance, Inspection, and Testing</b>	
<b>Reassembly Instructions</b>	
<b>Packaging</b>	
<b>Sterilization</b>	
<b>Storage</b>	
<b>Additional Information</b>	

<b>Related Healthmark Products</b>	
<b>Other Product Support Documents</b>	ProFormance™ Brochure, ProFormance™ Price List
<b>Reference Documents</b>	
<b>Customer Service Contact</b>	Healthmark Industries Company, Inc. 18600 Malyn Blvd. Fraser, MI 48026 1-586-774-7600 <a href="mailto:healthmark@hmark.com">healthmark@hmark.com</a> hmark.com