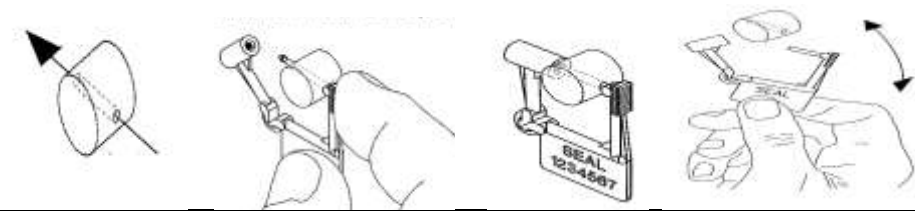


Brand Name of Product	Twist-Loks
Generic Name of Product	Tamper Evident Seals
Product Code Number(s)	5224
Purpose of Product	To show tampering and/or access to critical equipment and supplies.
Range of Applications for Product	Crash carts, emergency boxes, cabinets, and other areas where demonstrating tampering is critical.
Key Specifications of Product	<ul style="list-style-type: none"> • Numbered • Seal and broken by hand • Cannot be resealed once broken • Variety of colors available

Shipping & Storage

Shipping Conditions & Requirements	
Storage Conditions	
Packaging Contents	100 per box
Shelf Life	

Instructions for Using Product

Description of Use (s)	Seal crash carts, emergency boxes and other areas where evidence of tampering is critical
Preparation for Use	
Diagrams (drawings, pictures)	
Steps for Use of Product	<ol style="list-style-type: none"> 1. Make sure aperture is ready to accept seal. 2. Place arrow completely through aperture slot/opening. 3. Using thumb and fore finger, apply pressure and close. You should hear a click to ensure proper seal engagement. 4. Installed seal. 5. To remove, twist in any direction to break the shaft/chevron of seal.
Interpretation of Results	
Contraindications of Test Results	
Documentation	Follow the policy of your facility for recording date, staff and other information critical to management of critical inventory.
Special Warnings and Cautions	If a seal is broken, assume items have been tampered with or removed. Follow your facility instructions for inspection and replacement of critical inventory.
Disposal	

Reprocessing Instructions

Point of se	
Preparation for Decontamination	
Disassembly Instructions	
Cleaning – Manual	
Cleaning – Automated	
Disinfection	
Drying	
Maintenance, Inspection, and Testing	
Reassembly Instructions	
Packaging	
Sterilization	Compatible with steam sterilization. Standard cycle: 132°C (270°F) for 4 minutes
Storage	

Additional Information	
Related Healthmark Products	
Other Product Support Documents	Security Products Brochure, Security Products Price List
Reference Documents	
Customer Service Contact	Healthmark Industries Company, Inc 18600 Malyn Blvd. Fraser, MI 48026 1-586-774-7600 healthmark@hmark.com hmark.com

2020-11-24 msmith