

ProFormance™ Washer Monitoring Weekly Log Sheet: 1 Level Tunnel Washer (TWTK-1L)

Date:

Facility:

Washer: Rack No:

Name:

Water Quality
Aquachek
(circle result)

Target Values

p:

A:

H:

Cold Water



Comments

Hot Water




Comments

Pre-Wash Temperature
TempaChek-90

Temperature Should
Not Exceed 100°F

Record Result
(darken temperature reached)

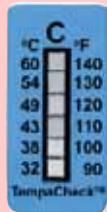

Circle Result

Pass

Fail

Comments

Verification of Cleaning Efficacy TOSI

Record Result
(circle result)



Comments

Thermal Disinfection
TempaChek-170

Target Temp:

Record Result
(circle highest temp achieved)



Comments

Cavitation Test Sono- Check Place one (1) test in sonic tank (circle resulting color)



Comments

	Bottom		Top		Sides		Bottom		Top		Sides		
	yes	no	yes	no	yes	no	yes	no	yes	no	yes	no	
Spray nozzles/arms are free of debris	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="text"/>
Nozzles(holes) properly aligned at target surface (up & down)	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="text"/>
All spray arms are present	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="text"/>
Spray arm spin freely	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="text"/>
Spray arm bushings are intact	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="text"/>
	yes	no			yes	no			yes	no			
Debris screen (in bottom of chamber) is clear of debris	<input type="checkbox"/>	<input type="checkbox"/>			<input type="checkbox"/>	<input type="checkbox"/>			<input type="checkbox"/>	<input type="checkbox"/>			<input type="text"/>
No staining/scaling from detergent, hardwater, etc.	<input type="checkbox"/>	<input type="checkbox"/>			<input type="checkbox"/>	<input type="checkbox"/>			<input type="checkbox"/>	<input type="checkbox"/>			<input type="text"/>
Detergent/enzyme at sufficient level in container	<input type="checkbox"/>	<input type="checkbox"/>											<input type="text"/>
Reviewed Cycle Settings (compared to Master)	<input type="checkbox"/>	<input type="checkbox"/>											<input type="text"/>
Other observations (comment)	<input type="text"/>											<input type="text"/>	
Recommended Actions	<input type="text"/>												
	<input type="text"/>												
	<input type="text"/>												