

ProFormance™ Washer Monitoring Log Sheet: 1 Level Tunnel Washer (TWTK-1L)

Date:

Facility:

Washer: Rack No:

Name:

Water Quality
Aquachek
(circle result)


Target Values

p:

A:


H:

Cold Water

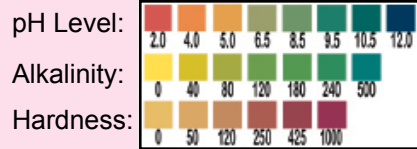

Comments

Hot Water

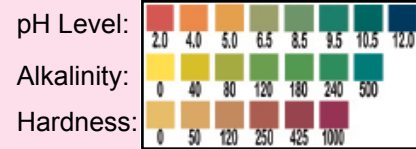

Comments

Pre-Wash Temperature
TempaChek-90

Temperature Should
Not Exceed 100°F



Record Result
(darken temperature reached)

Circle Result


Pass

Fail

Comments

**Verification
of
Cleaning
Efficacy
TOSI**

Record Result
(circle result)

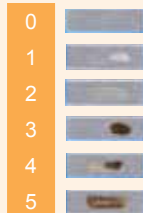


Comments

**Thermal
Disinfection**
TempaChek-170

Target Temp:

Record Result
(circle highest temp achieved)

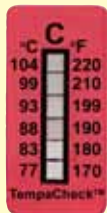

Comments

Cavitation Test
SonoCheck

Place one (1)
test in sonic tank
(circle resulting color)



Comments

Date:	Prewash Chamber			Wash Chamber			Disinfect Chamber			Comments			
	Bottom		Top		Sides		Bottom		Top		Sides		
	yes	no	yes	no	yes	no	yes	no	yes	no	yes	no	
Spray nozzles/arms are free of debris	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Nozzles(holes) properly aligned at target surface (up & down)	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
All spray arms are present	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Spray arm spin freely	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Spray arm bushings are intact	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
	yes no				yes no				yes no				
Debris screen (in bottom of chamber) is clear of debris	<input type="checkbox"/>	<input type="checkbox"/>			<input type="checkbox"/>	<input type="checkbox"/>			<input type="checkbox"/>	<input type="checkbox"/>			
No staining/scaling from detergent, hardwater, etc.	<input type="checkbox"/>	<input type="checkbox"/>			<input type="checkbox"/>	<input type="checkbox"/>			<input type="checkbox"/>	<input type="checkbox"/>			
Detergent/enzyme at sufficient level in container	<input type="checkbox"/>	<input type="checkbox"/>											
Reviewed Cycle Settings (compared to Master)	<input type="checkbox"/>	<input type="checkbox"/>											
Other observations (comment)	<input type="text"/>											<input type="text"/>	
Recommended Actions	<input type="text"/>												
	<input type="text"/>												
	<input type="text"/>												