

ProFormance™ Washer Monitoring Log Sheet: 3 Level Rack (AWTK-3L)

Date: Facility: Washer: Rack No: Name:

Water Quality
Aquachek
(circle result)

Target Values

p:

A:


H:

Cold Water



Comments


Hot Water

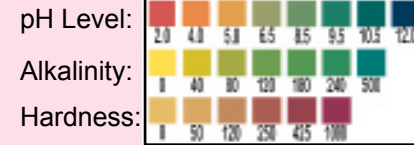

Comments

Pre-Wash Temperature
TempaChek-90

Temperature Should Not
Exceed 100°F

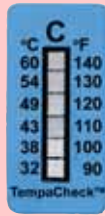

Record Result
(darken temperatures reached)

Circle Result

Pass

Fail

Comments

Verification
of
Cleaning
Efficacy
TOSI

Record Result
(circle result)


	Bottom	Middle	Top
0			
1			
2			
3			
4			
5			

Comments

Thermal Disinfection
TempaChek-170

Target Temp:

Record Result
(circle highest temp achieved)

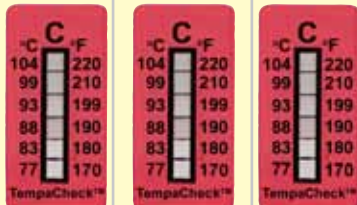

Comments

Date:	Chamber				Racks				Comments	
	Bottom		Top		Bottom		Top			
	yes	no	yes	no	yes	no	yes	no		
Spray nozzles/arms are free of debris	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>		
Nozzles(holes) properly aligned at target surface (up & down)	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>		
All spray arms are present	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>		
Spray arm spin freely	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>		
Spray arm bushings are intact	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>		
	yes	no								
Debris screen (in bottom of chamber) is clear of debris	<input type="checkbox"/>	<input type="checkbox"/>								
Instrument rack coupling with manifold properly	<input type="checkbox"/>	<input type="checkbox"/>								
No staining/scaling from detergent, hardwater, etc.	<input type="checkbox"/>	<input type="checkbox"/>								
Detergent/enzyme at sufficient level in container	<input type="checkbox"/>	<input type="checkbox"/>								
Reviewed Cycle Settings (compared to Master)	<input type="checkbox"/>	<input type="checkbox"/>								
Other observations (comment)	<input type="text"/>								<input type="text"/>	
Recommended Actions	<input type="text"/>									