

ProFormance™ Washer Monitoring Log Sheet: 4 Level Rack (AWTK-4L)

Date:

Facility:

Washer: Rack No:

Name:

Water Quality Aquachek (circle result)

Target Values

p:

A:

H:

pH Level:
Alkalinity:
Hardness:



Comments

Hot Water


pH Level:
Alkalinity:
Hardness:



Comments

Pre-Wash Temperature TempaChek-90

Temperature Should Not Exceed 100°F

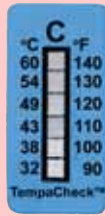

Record Result
(darken temperatures reached)

Circle Result

Pass

Fail

Comments

Verification of Cleaning Efficacy TOSI

Record Result
(circle result)

	Bottom	Middle L	Middle U	Top
0				
1				
2				
3				
4				
5				

Comments

Thermal Disinfection TempaChek-170

Target Temp:

Record Result
(circle highest temp achieved)



Comments

Date:	Chamber				Racks						Comments	
	Bottom		Top		Bottom		Middle		Top			
	yes	no	yes	no	yes	no	yes	no	yes	no		
Spray nozzles/arms are free of debris	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="text"/>
Nozzles(holes) properly aligned at target surface (up & down)	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="text"/>
All spray arms are present	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="text"/>
Spray arm spin freely	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="text"/>
Spray arm bushings are intact	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="text"/>
	<input type="checkbox"/>	<input type="checkbox"/>										
Debris screen (in bottom of chamber) is clear of debris	<input type="checkbox"/>	<input type="checkbox"/>										<input type="text"/>
Instrument rack coupling with manifold properly	<input type="checkbox"/>	<input type="checkbox"/>										<input type="text"/>
No staining/scaling from detergent, hardwater, etc.	<input type="checkbox"/>	<input type="checkbox"/>										<input type="text"/>
Detergent/enzyme at sufficient level in container	<input type="checkbox"/>	<input type="checkbox"/>										<input type="text"/>
Reviewed Cycle Settings (compared to Master)	<input type="checkbox"/>	<input type="checkbox"/>										<input type="text"/>
Other observations (comment)	<input type="text"/>											
Recommended Actions	<input type="text"/>											