


SteriScan™ Integrator parallels the performance of biological indicator spore strips with a safety margin throughout the entire sterilization range.

## Accurate and Instant Results with SteriScan.

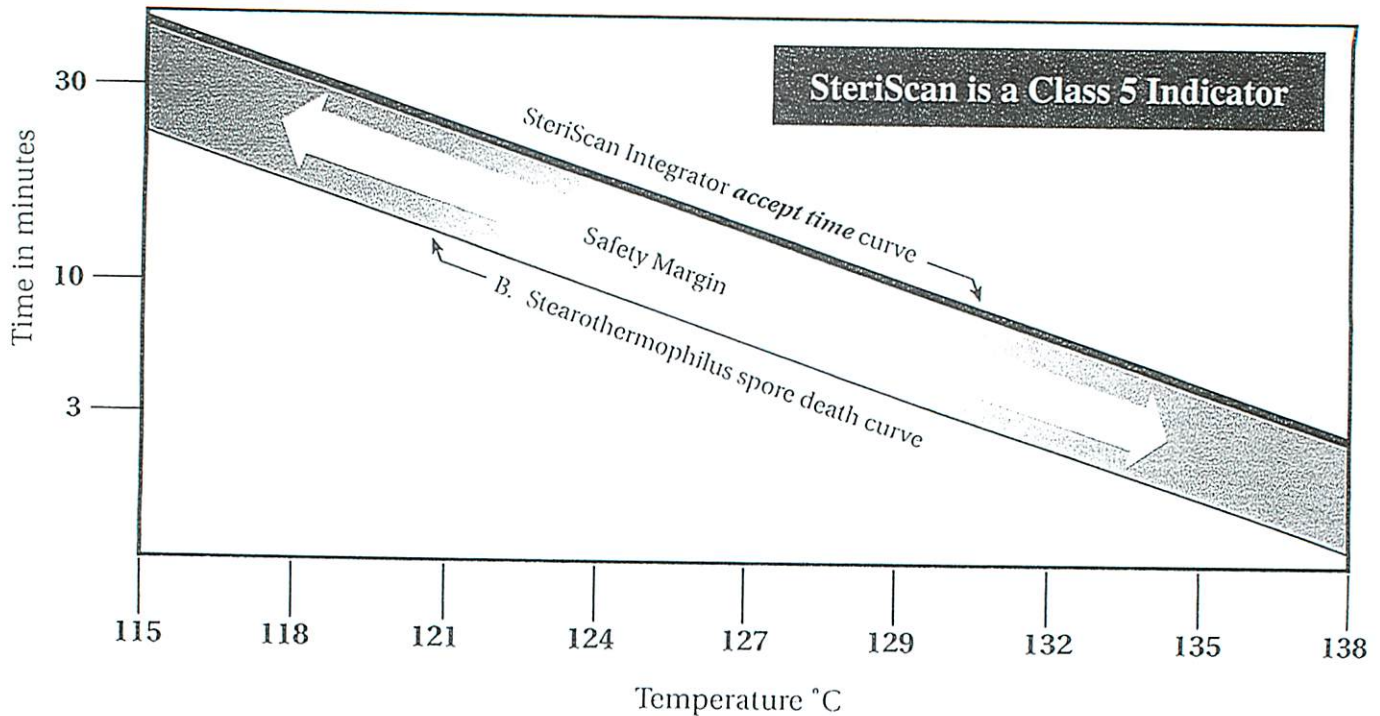

### Accurately Integrates All Critical Steam Sterilization Parameters

- Time
- Temperature
- Steam

### Quality You Can Depend On

- Developed by the inventor of the original integrator.
- Made with safe and non-toxic chemicals.

### Easy To Read

- Advancing color bar.
- Not a color changing strip.
- Eliminates the need to interpret or match colors.

### Effective Across The Normal Sterilization Range (115°C to 138°C)

- Eliminates the need to carry multiple inventories.

Distributed by:

DISTRIBUTED BY:  
HEALTHMARK INDUSTRIES, CO.  
33671 DOREKA DRIVE  
FRASER, MI 48026  
1-800-521-6224 / WWW.HMARK.COM

Reorder Number **SC1100**

Packaged 250 units in a foil bag and 8 bags (2,000 units) to a case.


# STERISCAN™

## STEAM STERILIZATION INTEGRATOR


### Takes The Guesswork Out Of Steam Sterilization




Unprocessed

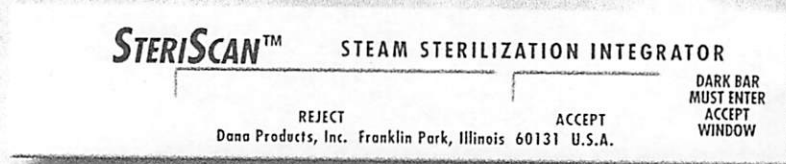

Processed  
Reject

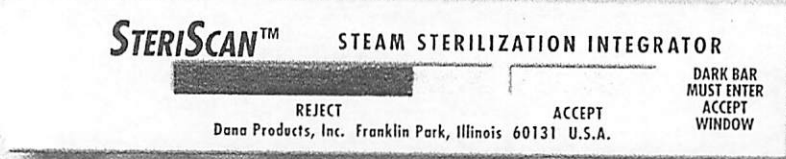

Processed  
Accept

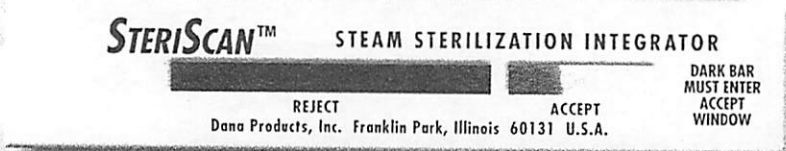

If the dark bar has not reached the *accept window* adequate sterilization conditions are not satisfied. Material should be reprocessed. *Product shown actual size: 3/4" x 4".*